

Parameters

Classification		Parameter	Measurement range			Measurement time (min.)	
			Unit (A)		Unit (B)		
Biochemical tests	Enzymes	ALP	14 – 1183	U/L	0.23 – 19.76	μ Kat/L	4
		v-AMY	100 – 2500	U/L	1.67 – 41.75	μ Kat/L	3
		CPK	10 – 2000	U/L	0.17 – 33.40	μ Kat/L	4
		GGT	10 – 1200	U/L	0.17 – 20.04	μ Kat/L	5
		GOT/AST	10 – 1000	U/L	0.17 – 16.70	μ Kat/L	4
		GPT/ALT	10 – 1000	U/L	0.17 – 16.70	μ Kat/L	4
		v-LIP	10 – 1000	U/L	0.17 – 16.70	μ Kat/L	5
		LDH	50 – 900	U/L	0.84 – 15.03	μ Kat/L	2
	General chemistry	ALB	1.0 – 6.0	g/dL	10 – 60	g/L	6
		BUN	5.0 – 140.0	mg/dL	1.79 – 49.98	mmol/L	4
		Ca	4.0 – 16.0	mg/dL	1.00 – 4.00	mmol/L	4
		CRE	0.2 – 24.0	mg/dL	18 – 2122	μ mol/L	5
		GLU	10 – 600	mg/dL	0.6 – 33.3	mmol/L	6
		IP	0.5 – 15.0	mg/dL	0.16 – 4.84	mmol/L	5
		Mg	0.2 – 7.0	mg/dL	0.08 – 2.88	mmol/L	4.5
		NH ₃	10 – 500	μg/dL	7 – 357	μ mol/L	2
		TBIL	0.2 – 30.0	mg/dL	3 – 513	μ mol/L	6
		TCHO	50 – 450	mg/dL	1.29 – 11.64	mmol/L	6
		TCO ₂	5 – 40	mmol/L	5 – 40	mmol/L	5
		TG	10 – 500	mg/dL	0.11 – 5.65	mmol/L	4
		TP	2.0 – 11.0	g/dL	20 – 110	g/L	6
		UA	0.5 – 18.0	mg/dL	30 – 1071	μ mol/L	4
	Electrolytes (on only one slide)	Na	75 – 250	mEq/L	75 – 250	mmol/L	1
		K	1.0 – 14.0	mEq/L	1.0 – 14.0	mmol/L	
		Cl	50 – 175	mEq/L	50 – 175	mmol/L	
Immunological Test		vc-CRP	0.3 – 7.0	mg/dL	3 – 70	mg/L	5
Comprehensive S-panel		TP, ALB, ALP, GLU, TBIL, IP, TCHO, GGT, GPT/ALT, Ca, CRE, BUN					
Pre-Surgical S-panel		TP, ALP, GLU, GPT/ALT, CRE, BUN					
Plus S-panel		TG, v-LIP, GOT/AST, v-AMY, Mg, Na-K-Cl					
Calculations		GLOB, ALB/GLOB, GOT/GPT, BUN/CRE, Na/K, Anion Gap					

* Unit (A) and (B) is available
* There are parameters which may not be available in your area. For details please contact your local distributor.
* v-LIP, v-AMY and vc-CRP slides are designed for veterinary use.
Other slides are designed for use in analyzing human blood. Measurement results may vary depending on the animal species. For details, please contact your local distributors.

Specifications

Measurement test	Colorimetry 23 tests Electrolytes 3 tests
Throughput	Colorimetry 120 test/ hour Combined 128 test/ hour
Number of sample rack	1
Number of incubator cell	Colorimetry 13, Electrolytes 1
Measurement time	Colorimetry 2 to 6 minutes/ test Electrolytes 1 minute/ 3 tests (Na, K, Cl)
Sample type	Plasma, Serum, Whole blood
Sample volume	Colorimetry 10μL/ test Electrolytes 50μL/ 3 tests (Na, K, Cl)
Data transmission to PC	RS 232C (1 port), USB (2 port), LAN (1 port)
Data print	Thermal printer
Electrical requirements	Single phase AC; 100 - 240 V ±10%; 50 to 60 Hz
Display	7 inch color touch panel
Dimensions	470 (W) 360 (D) 370 (H) mm
Weight	Approx. 27kg (NX600V IC) / Approx. 26kg (NX600V)
Operating temperature	15 to 32°C (59 to 89 °F)
Operating humidity	30 to 80% RH (no vapor condensation)
Altitude	Up to 2,000m (810 hPa)

Dimensions (mm)



NX600V Series Products

	With electrolyte measuring function	With automatic dilution function	With centrifugal separation function
NX600V IC	●	●	●
NX600V	●	●	—

Option Items

Barcode Reader
Barcode reader is available as option item to read sample ID on sample tube.



DRI-CHEM NX600V (Product name: FUJII DRI-CHEM NX600V IC, FUJII DRI-CHEM NX600V, FUJII DRI-CHEM NX600iV)

Read more about our Medical Solutions: www.fujifilm.com



MSO – Medical Solutions
Rua do Facho nº 9 Edif. Altamira Loja A 3770-056 Oitã - Portugal
+351 234 062 334 | +351 965 079 682 | geral@mso.pt

Ref. No. NX600 VET E FXX-XB-0400 (21•01•F1079•F8598) Printed in Japan ©2020 FUJIFILM Corporation

FUJIFILM
Value from Innovation

Small body, big flexibility.



10μL
/test

Automated Clinical Chemistry Analyzer **NEW**

DRI-CHEM NX600V
For veterinary

3-step easy measurement



Using the item code recorded on the slide, the item and lot number are automatically identified. After pressing the start key everything through the printing of the measurement data is fully automatic.


Test results
can be printed



One pack panels (S-panels)

Set multiple required items at once

One pack panels are also available. There is no need to open individually packed slides in a panel. This improves work efficiency and reduces waste.

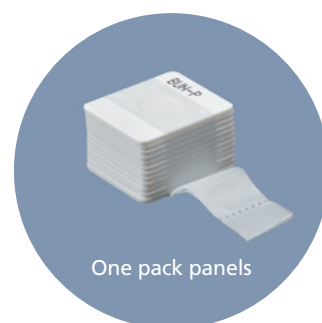

Individually packed slides in a panel



Simply place it in the slide cartridge and remove the binding tape



Measure directly from whole blood

Built-in centrifugal separation unit*

Simply insert the whole blood in the dedicated IC cup and place the cup in the unit. Automatically execute everything from plasma extraction to analysis. *NX600V IC



Can also measure serum and plasma separated by external centrifugation

This is possible when there is a need to reduce the amount of blood sampled or use the plasma sample in other tests.




Continuous measurement

Reduce waiting time

As soon as the pipetting procedure for the first sample is completed, the slides for the next sample can already be set and the test started. This minimizes waiting time for your patients.


NX600V processing cycle in minutes for 18 clinical parameters per patient sample

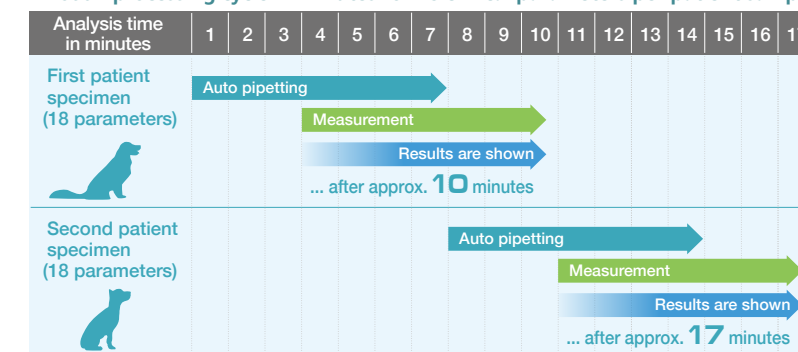

* In case of whole blood samples, it will take another 2 minutes for the centrifugation for each sample.

Improved design supports smooth workflow

User-friendly Fujifilm design

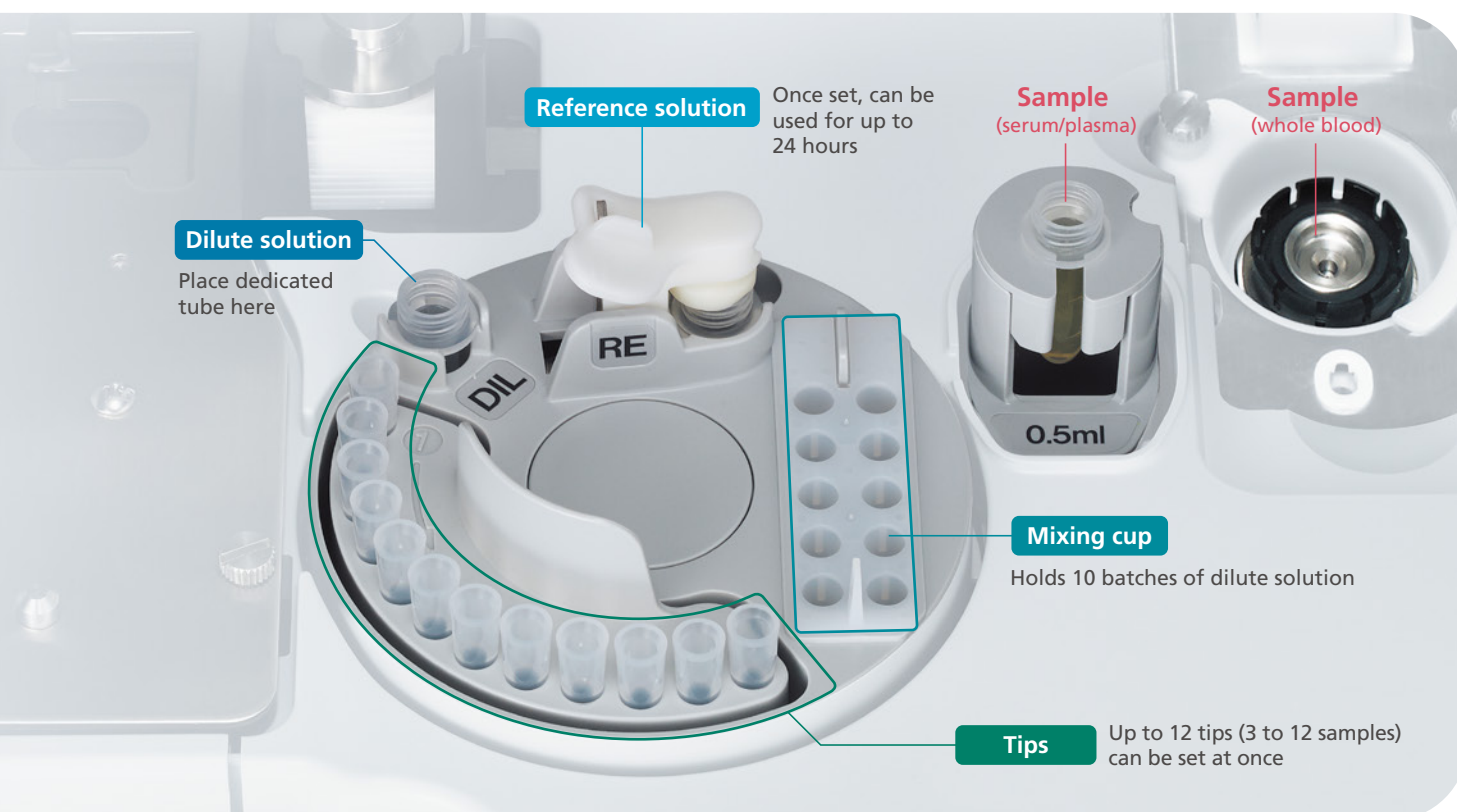


Set consumables for multiple samples at once

The layout is simpler than our previous models, and this makes it easier to set consumables. Multiple tips can be set at once, which eliminates the need to set reference solution for each measurement. This reduces the frequency of consumable replacement.



Large 7-inch intuitive LCD touch panel

The simple and easy-to-use operation screen provides a swift transition speed of 0.5 seconds or less for stress-free operation. A rich variety of help screens are also provided. This enables easy operation while checking illustrated usage and maintenance procedures.

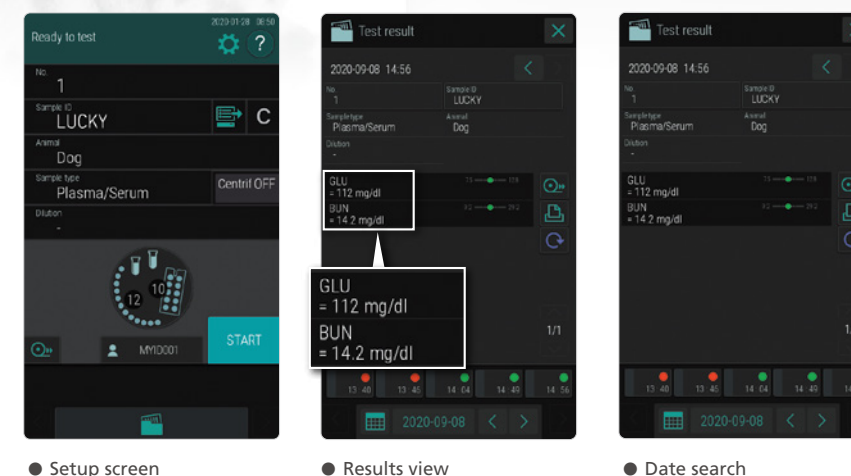

A rich variety of functions for flexible support



Microsample testing

Only 10 μ L of sample is required for each test. (CRP: 5 μ L/test, ISE: 50 μ L/test, dead volume is necessary using automatic measurement.) A very useful feature especially in small animals. If there is less sample, measurement can be done by manual pipetting.



High throughput

128 tests can be processed per hour.



Automatic dilution function

Labor intensive operations like dispensing, mixing etc. are automated. The only user operation required is when selecting the dilution ratio.



No Calibration

The FUJI DRI-CHEM system uses dry type reagents which are very stable. As the reagents are stable, calibration is not necessary and it utilizes a standard curve for the measurements. Variability in the reagent lots are adjusted using the QC card.

*ISE: Calibration and correction are not necessary. QC card is not attached.

Easy reagent lot information input by QC card

Information for every reagent lot that are needed for the measurement are easily input into the machine using the QC card.



Reduces biohazard risks

Used slides and tips fall automatically into a disposal box. Reduces the risk of contact with specimens.



Achieve high-precision testing with reliable technology

High-quality slides utilize coating technology cultivated in the field of photography.

FUJI DRI-CHEM SLIDE Two types of slides are provided for different test items.



Colorimetric method slide

Enzymes, General chemistry, and Immunology


This multilayered slide is composed of dry chemical ingredients needed for the reaction and other functional materials. It quantifies enzymes and chemicals using colorimetric method.




Potentiometric method slide

Electrolytes


Each slide comes with an ion selective film electrode for each of Na, K, and Cl. Slides quantify electrolytes in the sample by a potentiometric method.

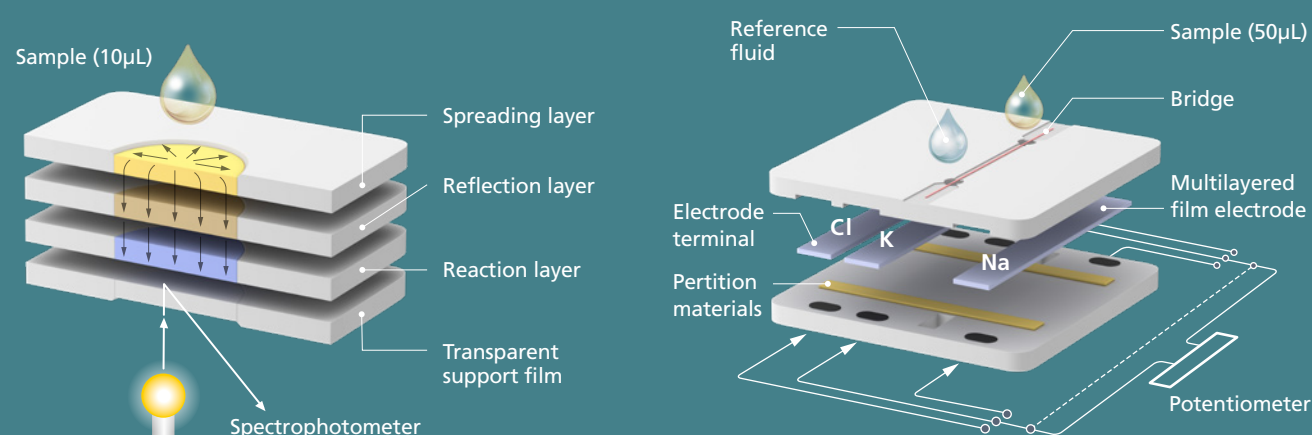

Centrifugal separation cup (IC cup)

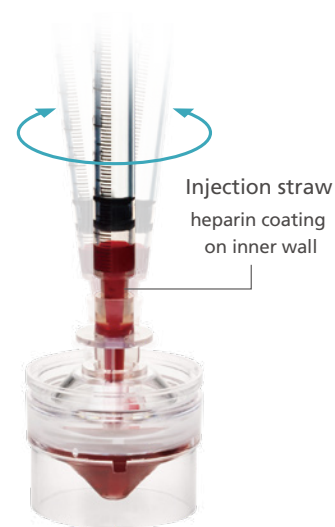


The straw can be easily removed simply by tilting it.

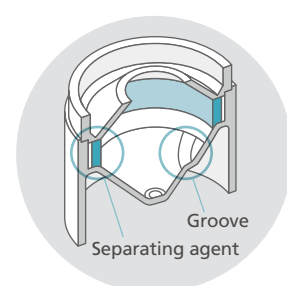


Unique Fujifilm structure reduces risks

Inserting a syringe into the straw included with the cup and injecting the sample prevents sample scattering. Furthermore, an internal structure with a separation layer and groove reduces hemolysis and promotes the stable supply of good quality plasma.

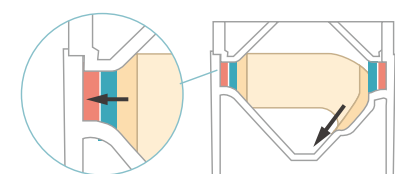

Centrifugal separation traps the blood cells with a large specific gravity on the outer side of the separating agent.


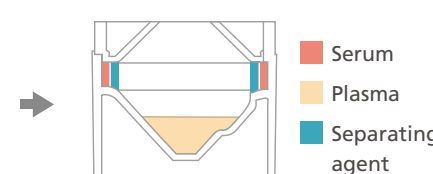


After separation, the plasma follows the groove and accumulates separately in the bottom of the cup.

■ Serum
■ Plasma
■ Separating agent


Based on technology and experience developed over many years, Fujifilm DRI-CHEM provides conclusive test results

Fujifilm DRI-CHEM slides achieve an extremely high degree of precision thanks to a special filter layer that reliably filters out any substances that could interfere with the serum or plasma specimen. The high quality of Fujifilm reagents is demonstrated by a comparison with the standard measurement methods.

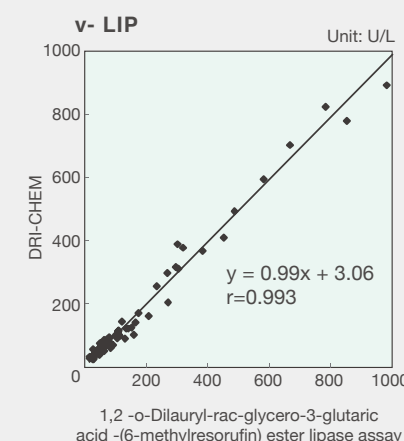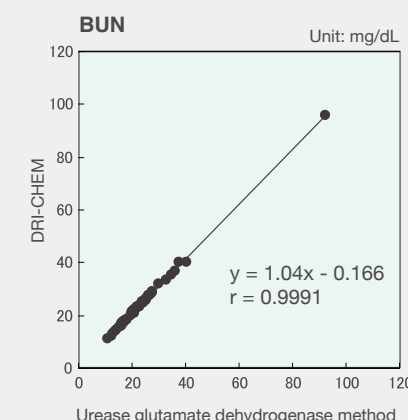


Realtime diagnosis canine CRP test for inflammation screening

NX700V can measure lipase activity for canine CRP with the rest of the biochemistry test items.

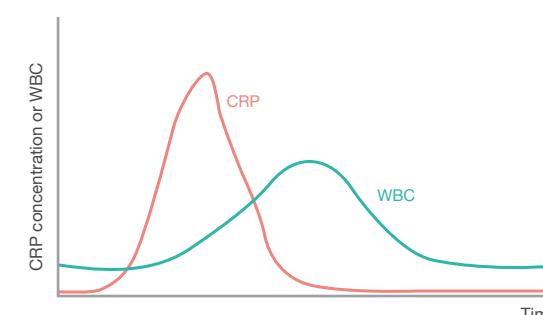

NX600V can measure lipase activity for canines

Lipase activity, which is useful for the screening of pancreatitis, can also be easily measured. This contributes to quick screening of pancreatitis onsite.

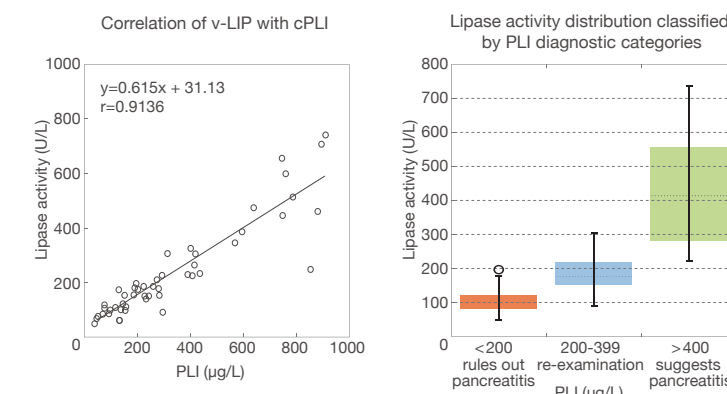

Citation: J. Vet. Med. Sci. 73(11):1481-1483, 2011 Katsumi Ishioka
J-vet vol.291 2011.6 pp13-15, pp25-29

DRI-CHEM product lineup providing total in-house diagnostic testing



DRI-CHEM NX700V
Clinical chemistry analyzer

• High throughput 190 test/ hour (colorimetry and ISE combined)
• 5 sample rack • Without centrifugal separation unit



DRI-CHEM NX600V
Clinical chemistry analyzer

• Single sample rack



DRI-CHEM IMMUNO AU10V
Immunoassay analyzer

v-T4, vc-TSH, v-COR, v-BA, v-PRG, vf-SAA